

SmartRedux

SEPAX™ APPLICATION SOFTWARE

Our application software provides a fully-automated method to separate cellular product into its components in combination with our Sepax™ Cell Separation System. It enables a reproducible and user independent procedure.

The SmartRedux application is intended for concentration of the total nucleated cell (TNC) from bone marrow.


Key features

| Feature | Description |
|-------------------------|---------------|
| Initial volume | 30 to 3300 mL |
| Selectable final volume | up to 1500 mL |
| Average TNC recovery | 75%* |

*Indication only. Performance depends on user/cellular product configuration.

Workflow

| Step | Details |
|------|--|
| 1 | After bone marrow aspiration and filtration, connect initial product to the kit. |
| 2 | Install the kit on the Sepax™ instrument. |
| 3 | The Sepax™ system automatically performs the procedure. |
| 4 | Collect the concentrated cells from the final product bag. |

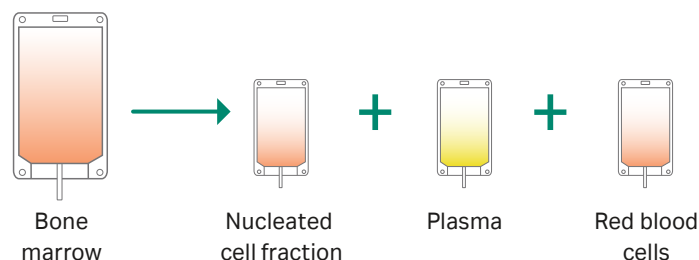

Product specifications

- Fixed or proportional option available for the final volume, to fit different needs
- Additional plasma parameter to optimize final nucleated cell fraction characteristics
- Traceability IDs: Use the Sepax™ traceability kit (barcode reader + printer) to track specific identification criteria and print out a PDF report
- Main user interface available in multiple languages

Ordering Information

| Products | Product code |
|---------------------------------------|--------------|
| SmartRedux Sepax™ 2 protocol software | 14301 |

Packaging: One box containing a USB stick with the installer.

| Related products | Product code |
|--|--------------|
| Sepax™ 2 S-100 | 14000 |
| Sepax™ 2 RM | 14100 |
| CS-490.1 Sepax™ cell separation kit - Box of 6 Units | 10001 |

Not all products are available in all countries. Please contact your sales representative about the product availability in your region.

Manufacturer information

Biosafe SA, Route du Petit-Eysins 1, 1262 Eysins, Switzerland

Contact us at:

www.cytivalifesciences.com/en/us/support/contact-us

cytiva.com

Cytiva and the Drop logo are trademarks of Global Life Sciences IP Holdco LLC or an affiliate. Sepax is a trademark of Global Life Sciences Solutions USA LLC or an affiliate doing business as Cytiva.

© 2021 Cytiva

All goods and services are sold subject to the terms and conditions of sale of the supplying company operating within the Cytiva business. A copy of those terms and conditions is available on request. Contact your local Cytiva representative for the most current information.

For local office contact information, visit cytiva.com/contact

CY21160-08Sep21-DF

